# ASEPTIQUIK® STC SERIES CONNECTOR



AseptiQuik® STC Connectors integrate the AseptiQuik® sterile connector and the Steam-Thru® II SIP connector, giving manufacturers even greater flexibility for hybrid stainless steel and single-use processing. The AseptiQuik STC connector combines either an AseptiQuik G or AseptiQuik C sterile connector with a Steam-Thru II connection that can be mounted directly to a stainless steel vessel via a sanitary termination. AseptiQuik STC connectors utilizing an AseptiQuik G sterile connector offer a genderless design, simplifying system integration at the aseptic connector end. AseptiQuik STC connectors utilizing an AseptiQuik C sterile connector allow system designers to incorporate a keyed female-to-male connection at the aseptic connector end.

The union of the two connectors into a single unit through a sanitary clamp allows an AseptiQuik sterile connection to be steamed on to stainless equipment via SIP. After the SIP cycle, a wide range of single-use systems can be connected. The SIP process can be done in advance allowing a quick and easy sterile connection to the AseptiQuik half without having to wait 30-60 minutes for SIP prior to harvest.

FEATURES	BENEFITS
AseptiQuik design	Innovative three-step connection process reduces risk of operator error
Genderless design	Eases integration of single-use systems with universal mating between AseptiQuik G Connectors
Innovative three-port steam design	Allows a true steam-through SIP process which eliminates "dead legs"
Robust construction	Repeatable and reliable performance with no additional hardware required
Patented steam valve design	Allows sterile connection and disconnection to stainless equipment and permits a high media flow rate.
Sanitary interface between the two connectors	More secure connection than tubing with cable ties
CPC Click	Audible confirmation of assembly steps
Market availability	Open access through multiple supply chain partners

## Specifications • • •

### 113

### PRESSURE:

**Steam position:** Up to 35 psi, 2.4 bar **Flow position:** Up to 20 psi, 1.4 bar

### TEMPERATURE:

#### Steam position:

Up to 275°F (135°C) for 60 minutes **Flow position:** 39°F to 104°F (4°C to 40°C)

#### **STERILIZATION:**

Gamma: Up to 50kGy irradiation

AutoClave: High Temp (HT) version: Up to 266°F (130°C) for 30 minutes (AQSTC) Up to 266°F (130°C) for 60 minutes (AQGSTC)

**Note**: A slight clockwise rotation of the clamp nut may be needed after autoclave.

#### **MATERIALS:**

#### **Main Components:**

AseptiQuik - Polycarbonate (white), USP Class VI, ADCF

Steam-Thru II - Polysulfone (amber tint), USP Class VI. ADCF

**Seals:** Silicone (clear), platinum-cured, USP Class VI, ADCF

#### **Pull Tahs**

Polycarbonate (blue, standard version), USP Class VI, ADCF Polycarbonate (white, HT version), USP Class VI, ADCF

Caps: Polypropylene (clear), USP Class VI, ADCF

#### Memhran


Polyethylene (standard version), USP Class VI, ADCF

Hydrophobic polyethersulfone (HT versions), USP Class VI, ADCF, PTFE strip sticker

Clamp: Nylon 66 (white), USP Class VI



# ASEPTIQUIK® STC SERIES DIMENSIONS

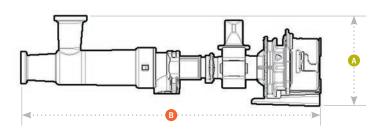

A = Height/Diameter
B = Total Length

• = Actuated Length

Mates with gendered male and female AseptiQuik Series.

## **Coupling Bodies**

• POLYCARBONATE with Blue Pull Tabs



PART NO. AQSTC2330900	A 2.78 (70.6 mm)	B 9.25 (235 mm)	8.42 (213.9 mm)	
AQSTC2331000	2.78 (70.6 mm)	9.25 (235 mm)	8.42 (213.9 mm)	

• POLYCARBONATE HT with White Pull Tabs




3/4" X 3/4" SANITARY STEAM-THRU II WITH ASEPTIQUIK BODY (PICTURED)

3/4" X 1-1/2" SANITARY STEAM-THRU II WITH ASEPTIQUIK BODY

ART NO. QSTC2330900HT	A 2.78 (70.6 mm)	B 9.25 (235 mm)	© 8.42 (213.9 mm)
QSTC2331000HT	2.78 (70.6 mm)	9.25 (235 mm)	8.42 (213.9 mm)

# ASEPTIQUIK® G STC SERIES DIMENSIONS

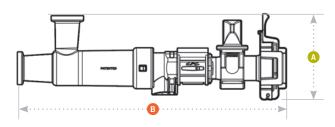

A = Height/DiameterB = Total LengthC = Actuated Length

## **Coupling Bodies**

• POLYCARBONATE with Blue Pull Tabs



TERMINATION
3/4" X 3/4" SANITARY STEAM-THRU II
WITH GENDERLESS ASEPTIQUIK G BODY

WITH GENDERLESS ASEPTIQUIK G BODY
3/4" X 1-1/2" SANITARY STEAM-THRU II
WITH GENDERLESS ASEPTIQUIK G BODY

PART NO. AQGSTC2330900	3.05	B 8.26	7.36
AQGSTC2331000	3.05	8.26	7.36

Mates with

AseptiQuik G

POLYCARBONATE HT with White Pull Tabs



TERMINATION
3/4" X 3/4" SANITARY STEAM-THRU II
WITH GENDERLESS ASEPTIQUIK G BODY

3/4" X 1-1/2" SANITARY STEAM-THRU II WITH GENDERLESS ASEPTIQUIK G BODY

PART NO. AQGSTC2330900HT	3.05	B 8.26	7.36
AQGSTC2331000HT	3.05	8.26	7.36

Smart fluid handling to take you forward, faster. 133